

Ace Glass Introduces...




A New Concept in Lab Filtration Apparatus!

**(All Borosilicate Glass
& Nylon Components)**


No clamp, no leakage!

Using Ace-Thred design, filter funnel and filter flask thread together without the silicone stopper and troublesome clamping device.

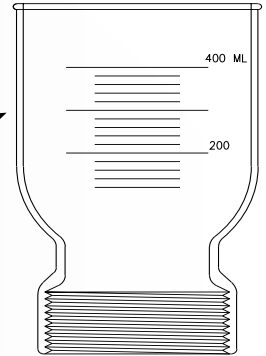



Easy cleaning, removable frits!

Unit disassembles completely for easy cleaning. Uses removable frits of different porosities, from extra coarse to extra fine, instead of expensive, sealed-in frits.

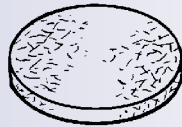

Better membrane utilization!

Accommodates the 47mm and 25mm membranes and uses 45mm of the 47mm grid vs. only 35-40mm in older designs.

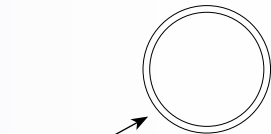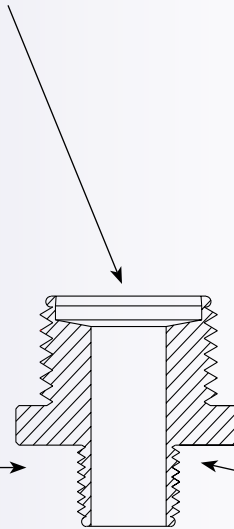

Nylon membrane retaining ring completes seal.

Choose your adapter!

Features autoclaveable nylon adapters, but you can also special order versions in PTFE or other polymers.



3702-10

Complete Filtration Apparatus — 47mm ♠

Description	Size	Order Code	Each
Includes: 1000mL Filter Flask w/Ace-Thred top, & #7 Ace-Thred & Ace-Safe connector, Nylon Adapter w/#25 Ace-Thred bottom, & #50 Ace-Thred top, nylon ring 25mm Porosity B (70-100 micron) fritted disc, 500mL graduated funnel w/#50 Ace-Thred bottom.	47mm	3702-10	353.78
1000mL Filter Flask w/#25 Ace-Thred & #7 Ace-Safe connector	1000mL	3702-08	127.68
1000mL Filter Flask w/#24/40 ST outer joint top and hose barb side port	1000mL	6979-15	107.41
Nylon adapter w/#25 Ace-Thred bottom & #50 Ace-Thred top	Adapter	3702-06	131.13
Nylon adapter w/#50 Ace-Thred top and 24/40 inner joint bottom	Adapter	3702-05	121.87
500mL Funnel w/#50 Ace-Thred bottom, graduated	500mL	3702-04	56.89
Porosity B fritted disk (70-100 micron) 47mm	47mm	3703-47	27.57
Porosity A fritted disk (145-174 micron) 47mm	47mm	3703-45	27.57
Porosity C fritted disk (25-50 micron) 47mm	47mm	3703-49	27.57
Nylon Ring-47mm	47mm	3702-02	10.51

♠ Subject to Ace Discount